

Evidence-Bench™

Designed for Easy Evidence Processing

Evidence-bench allows for mobile forensic evidence processing in your lab, garage or evidence room.

72



— Model EVB-72

Shown with optional airflow pod, magnifying lamp and power package outlet strip.

INTRODUCTION

The Air Science™ Evidence-Bench series of laboratory work tables are designed specifically for high volume processing of forensic evidence by multiple users. Professional design and construction features offer convenience during use, cleaning, and maintenance.

STANDARD FEATURES

- Kraft paper dispenser
- Lockable storage compartments to secure tools
- Open storage shelf
- Wheels allowing the bench to be moved throughout the building wherever processing needs to take place.
- Optional airflow pod

AIRFLOW POD

The optional drop-in Airflow Pod permits efficient processing of evidence when using powders and/or chemicals.

- Based on the Air Science Multiplex Filtration technology, the airflow pod is a high efficiency containment system that protects the user and the environment from


hazardous vapors or powders generated at the stainless steel work surface with down flow air movement.

- Unrestricted front and side access facilitates applications requiring complex and intensive involvement.
- The Airflow Pod eliminates putrid odors that may be emitted from the evidence being processed.



Made in the U.S.A.



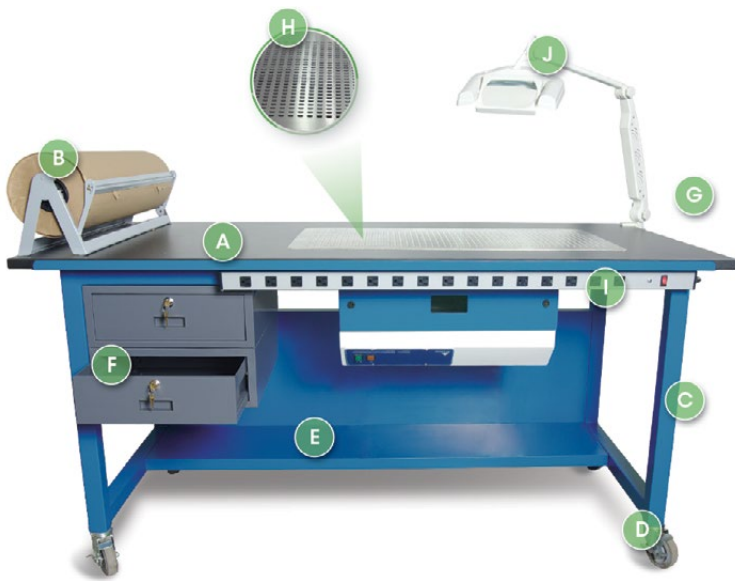

PRODUCT FEATURES:


- A. Epoxy resin work top chemical resistant
- B. Evidence kraft paper dispenser
- C. Heavy duty frame
- D. Locking wheels, heavy duty, non-marking
- E. Lower open shelf for storage
- F. Two locking supply compartments
- G. Open access from multiple sides
- H. Optional Airflow Pod filtration shown; see Accessories
- I. Optional 16 outlet power strip shown; see Accessories
- J. Optional magnifying lamp shown; see Accessories

• EVB-72 shown with accessories.

MODEL	DIMENSIONS		WEIGHT (lbs/Kg)	
	External (W x D x H)	Shipping (W x D x H)	Net	Ship
Evidence-Bench™				
EVB-72	72" x 28" x 36" 1828 x 711 x 914 mm	80" x 40" x 45" 2032 x 1016 x 1143 mm	325 / 147	400 / 182

PRODUCT SPECIFICATIONS

Evidence-Bench™ Model	EVB-72
Construction	Epoxy coated steel frame
Locking Wheels	Non-marking heavy duty caster wheels
Airflow (optional)	135.9 cfm
Work Surface	Black epoxy resin, chemical resistant
Kraft Paper Dispenser	Standard
Shelf	Standard
Drawers	Two, key locking, left hand side mount standard

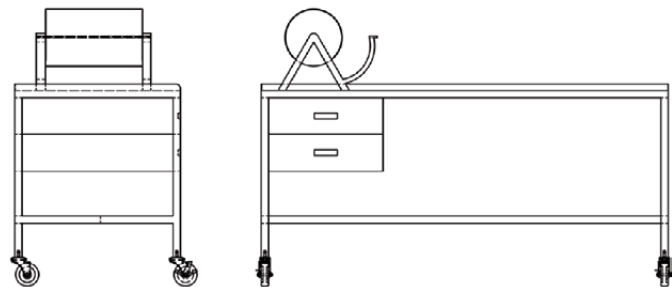

Side View

Evidence-Bench™ EVB-72

STANDARDS AND COMPLIANCE

Quality Management Systems	ISO 9001:2008
Environment	ISO 14001:2004 ENERGY STAR® Partner
Cabinet Performance	IEST-RP-CC002.2 AS 2252.6
Air Quality	ISO 14644-1, Class 4
Filtration	IEST-RP-CC034.1 IEST-RP-CC001.3 IEST-RP-CC007.1 EN 1822
Electrical Safety	UL-C-61010-1 CE Mark RoHS Exempt under EEE Category 9

OPTIONS AND ACCESSORIES

Drop-in Airflow Pod	P5-36H
Magnifying Lamp	MAG-L
Power Package Outlet Strip	PSTRIP-16, 16 outlets
Water Package	SINK-P, includes cup sink and water faucet
Stainless steel work surface	SS-TOP
Adjustable Levelling feet	LEVEL

* Factory installed. Specify when ordering.
Specifications are subject to change without notice.
Contact Air Science for custom sizes and options not listed herein.



PO Box 62296 • Fort Myers, FL 33907
T/239.489.0024 • Toll Free /800.306.0656 • F/800.306.0677
www.airscience.com

